



محصولات مراقبت قلبی Teleflex®

تهران، خیابان شهید قندی (پالیزی)، شماره ۱۲۸
کد پستی: ۱۵۵۴۷۴۳۷۶، تلفن: ۸۳۹۶ (خط ۳۰)
No.128, Qandi (Palizi) St., Tehran, Iran
Postal Code:1554744371, Tel: (+9821) 8396
www.iranpharmis.org



Superb Crossability in Challenging Cases

این خانواده از میکروکاتترها با هدف رساندن ابزار Intervention به تصلب های نرم و غیرآهکی در زوایای متنوعی طراحی شده است.


میکروکاتترهای SuperCross به علت داشتن Tip های زاویه دار امکان دسترسی به تصلب های عروق پرپیچ و خم و منشعب را برای پزشک فراهم می کنند.

شفاف اصلی با کوئل دولایه جهت تامین قدرت چرخش و پیچش بالا، انعطاف پذیری، پیش روندگی و مقاومت در برابر Kink شدن


لایه داخلی با پوشش پلیمر PTFE جهت تسهیل حرکت آسان و روان ابزار Intervention

پوشش هیدروفلیک در قسمت Distal میکروکاتتر جهت افزایش قدرت پیش روندگی و سهولت حرکت در مسیرهای پرپیچ و خم عروقی


غلاف پلاتینیوم/تنگستن در Tip جهت ردیابی Tip در نزدیکی تصلب

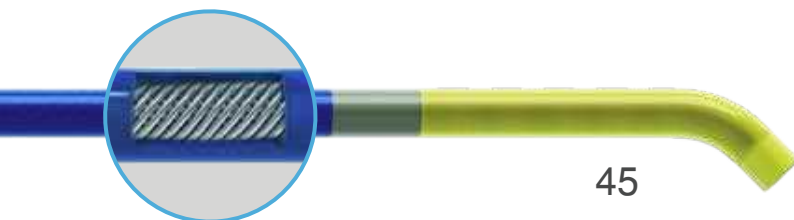

45

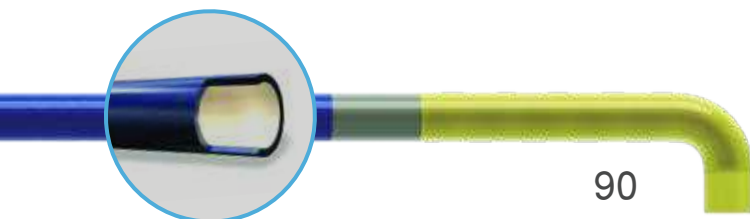

90

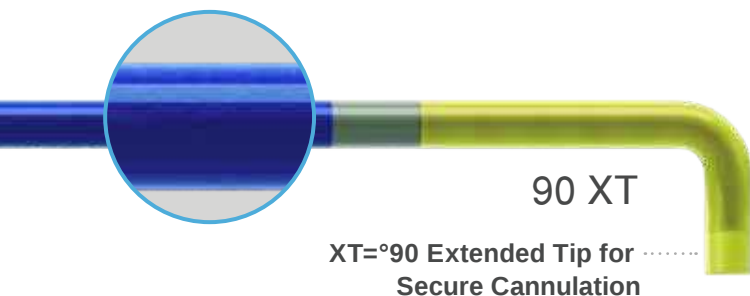

90 XT

XT=°90 Extended Tip for
Secure Cannulation



120

ویژگی ها:

- مناسب برای کیس های Retrograde و Antegrade
- مناسب برای عروق پریفرال و کرونری
- دارای قدرت پیش روندگی بالا به علت داشتن پوشش پلیمری یکنواخت و هیدروفیل
- دارای شفاف چند لایه سبک و انعطاف پذیر جهت دسترسی راحت تر به تصلب های غیرآهکی
- مناسب جهت دسترسی به عروق پر پیچ و غیر مستقیم به علت Tip زاویه دار


Ordering Information

SuperCross™ Microcatheters

The SuperCross™ Microcatheter is intended to be used in conjunction with steerable guidewires to access discrete regions of the coronary and/or peripheral vasculature. It may be used to facilitate placement and exchange of guidewires and other interventional devices and to subselectively infuse/deliver diagnostic and therapeutic agents.

MODEL	GUIDEWIRE COMPATIBILITY	DISTAL TIP	WORKING LENGTH	SHAFT CONSTRUCTION	DISTAL O.D.	PROXIMAL O.D.	DISTAL I.D.	PROXIMAL I.D.	HYDROPHILIC COATING
5302	0.014"	45	130 cm	Coiled	2.4 Fr. 0.79 mm / 0.031"	3.2 Fr. 1.07 mm / 0.042"	0.017" (0.43 mm)	0.018" (0.46 mm)	Distal 80 cm
5303	0.014"	45	150 cm	Coiled	2.4 Fr. 0.79 mm / 0.031"	3.2 Fr. 1.07 mm / 0.042"	0.017" (0.43 mm)	0.018" (0.46 mm)	Distal 80 cm
5304	0.014"	90	130 cm	Coiled	2.4 Fr. 0.79 mm / 0.031"	3.2 Fr. 1.07 mm / 0.042"	0.017" (0.43 mm)	0.018" (0.46 mm)	Distal 80 cm
5305	0.014"	90	150 cm	Coiled	2.4 Fr. 0.79 mm / 0.031"	3.2 Fr. 1.07 mm / 0.042"	0.017" (0.43 mm)	0.018" (0.46 mm)	Distal 80 cm
5306	0.014"	120	130 cm	Coiled	2.4 Fr. 0.79 mm / 0.031"	3.2 Fr. 1.07 mm / 0.042"	0.017" (0.43 mm)	0.018" (0.46 mm)	Distal 80 cm
5307	0.014"	120	150 cm	Coiled	2.4 Fr. 0.79 mm / 0.031"	3.2 Fr. 1.07 mm / 0.042"	0.017" (0.43 mm)	0.018" (0.46 mm)	Distal 80 cm
5308	0.014"	90 XT	130 cm	Coiled	2.4 Fr. 0.79 mm / 0.031"	3.2 Fr. 1.07 mm / 0.042"	0.017" (0.43 mm)	0.018" (0.46 mm)	Distal 80 cm
5309	0.014"	90 XT	150 cm	Coiled	2.4 Fr. 0.79 mm / 0.031"	3.2 Fr. 1.07 mm / 0.042"	0.017" (0.43 mm)	0.018" (0.46 mm)	Distal 80 cm

Packaged in quantities of 1 unit per box.

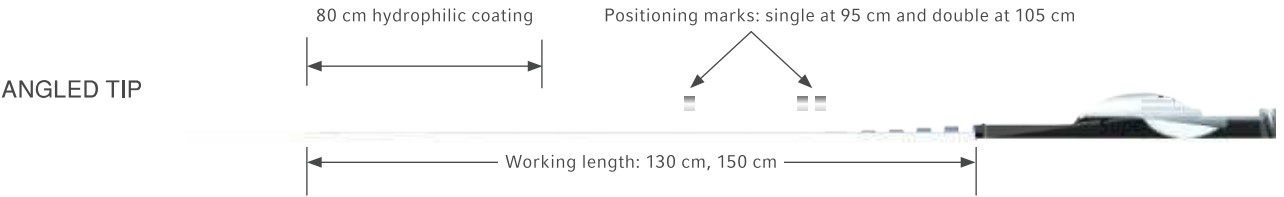

SuperCross™ Microcatheters—Straight| FT

The SuperCross™ Microcatheter is intended to be used in conjunction with steerable guidewires to access discrete regions of the coronary and/or peripheral vasculature. It may be used to facilitate placement and exchange of guidewires and other interventional devices and to subselectively infuse/deliver diagnostic and therapeutic agents.

MODEL	GUIDEWIRE COMPATIBILITY	DISTAL TIP	WORKING LENGTH	SHAFT CONSTRUCTION	DISTAL O.D.	PROXIMAL O.D.	DISTAL I.D.	PROXIMAL I.D.	HYDROPHILIC COATING
5300	0.014"	Straight	130 cm	Braided	1.8 Fr. 0.61 mm / 0.024"	2.5 Fr. 0.84 mm / 0.033"	0.017" (0.43 mm)	0.021" (0.53 mm)	Distal 40 cm
5301	0.014"	Straight	150 cm	Braided	1.8 Fr. 0.61 mm / 0.024"	2.5 Fr. 0.84 mm / 0.033"	0.017" (0.43 mm)	0.021" (0.53 mm)	Distal 40 cm
5340	0.014"	FT	130 cm	Braided	1.8 Fr. 0.61 mm / 0.024"	2.5 Fr. 0.84 mm / 0.033"	0.017" (0.43 mm)	0.021" (0.53 mm)	Distal 40 cm
5341	0.014"	FT	150 cm	Braided	1.8 Fr. 0.61 mm / 0.024"	2.5 Fr. 0.84 mm / 0.033"	0.017" (0.43 mm)	0.021" (0.53 mm)	Distal 40 cm

Packaged in quantities of 1 unit per box.
FT= Flexible Tip

